




VETPUMP® 2

Precision Suction & Irrigation System

Cleaning ears has never been so easy

The VETPUMP® 2 from KARL STORZ puts the safest and most effective flushing and suction at your fingertips!

More powerful, more precise and safer than other units, the VETPUMP® 2 is the ultimate device for lavage and suction during otoscopy, gastroscopy and more.

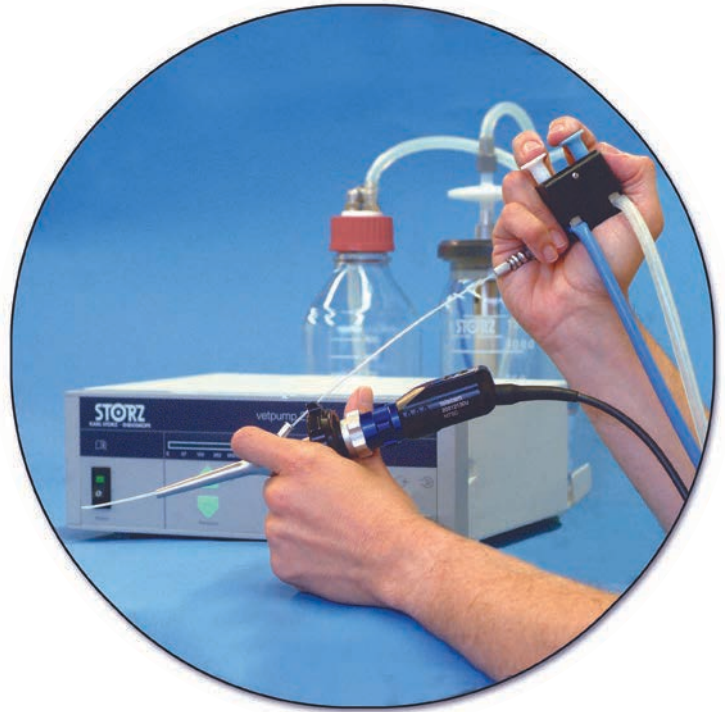

Features:

- Powerful suction and irrigation pump for both rigid and flexible KARL STORZ endoscopes
- There is no exchange of air between bottles, eliminating any potential contamination of clean irrigation fluids
- Micro-processor controlled unit is quiet and responsive; unit is completely silent when trumpet valves are not activated
- Irrigation pressure and suction vacuum can be regulated separately
- User-selected pressures are precise, powerful and reproducible
- Versatile hand piece attaches to catheters, irrigation tubes, or directly to endoscope channel
- Bottles, tubing and handle are fully autoclavable



69 3216 01

69 3216 01

VETPUMP® 2, power supply 100/120/230/240 VAC, 50/60 Hz including:



Suction Bottle, 1.5 l, sterilizable
Bottle Cap, for use with suction bottle 1.5 l, sterilizable




Irrigation Bottle, 1 l, sterilizable
Bottle Cap, for use with irrigation bottle

Bottle Holder, for suction bottle
Bottle Holder, for irrigation bottle
Bacterial Filter

Mains Cord
Mains Cord, US-version

Pressure Regulated (mmHg)	Flow Regulated (l/ min)	Suction Pressure (mmHg)	Power Supply	Dimensions (w x h x d) mm	Weight (kg)	Certified to
max. 600	max. 3.5	regulated: max. 580	100/120/230/240 VAC, 50/60 Hz	305 x 110 x 260	3.8	EN 61010-1

Accessories for Otoscopy:

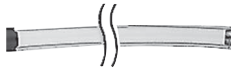

69 3200 33 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**




69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle



69810 **Handle**, for irrigation and suction, autoclavable, for use with tubing set 69811, in conjunction with 5 mm suction and irrigation tubes or catheters any size with LUER Lock connection



Accessories for Abdominal/Thoracic Lavage:

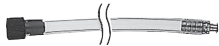

69 3200 33 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**



69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle




37112 A **Handle**, straight, with clamping valve, for suction and irrigation, **autoclavable**, for use with 5 and 10 mm suction and irrigation tubes



69811 LAP **Adaptor**, for use with suction and irrigation Handle 37112 A, plus Instrument Tubing Set 69811, package of 5



031121-10* **Gas filter**, for single use, sterile, package of 10, with male LUER-Lock connector for attachment to irrigation connector on unit

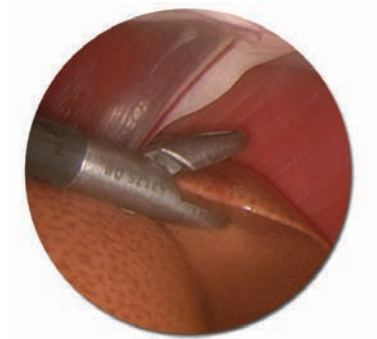

Accessories for Gastrointestinal Endoscopy:



20 3000 42 **Silicone Tubing Set**, for suction only, sterilizable



610 VE **Silicone Tubing Set**, for insufflation, attaches to light source adaptors 487 VE or 69487 KVE




*



mtp medical technical promotion gmbh, take-off GewerbePark 46, 78579 Neuhausen ob Eck/Germany, Tel.: +49 (0)7467 94504-0, Fax: +49 (0)7467 94504-99, E-Mail: info@mtp-tut.com, www.mtp-tut.com

VETPUMP® 2

System Components for use in Otoscopy

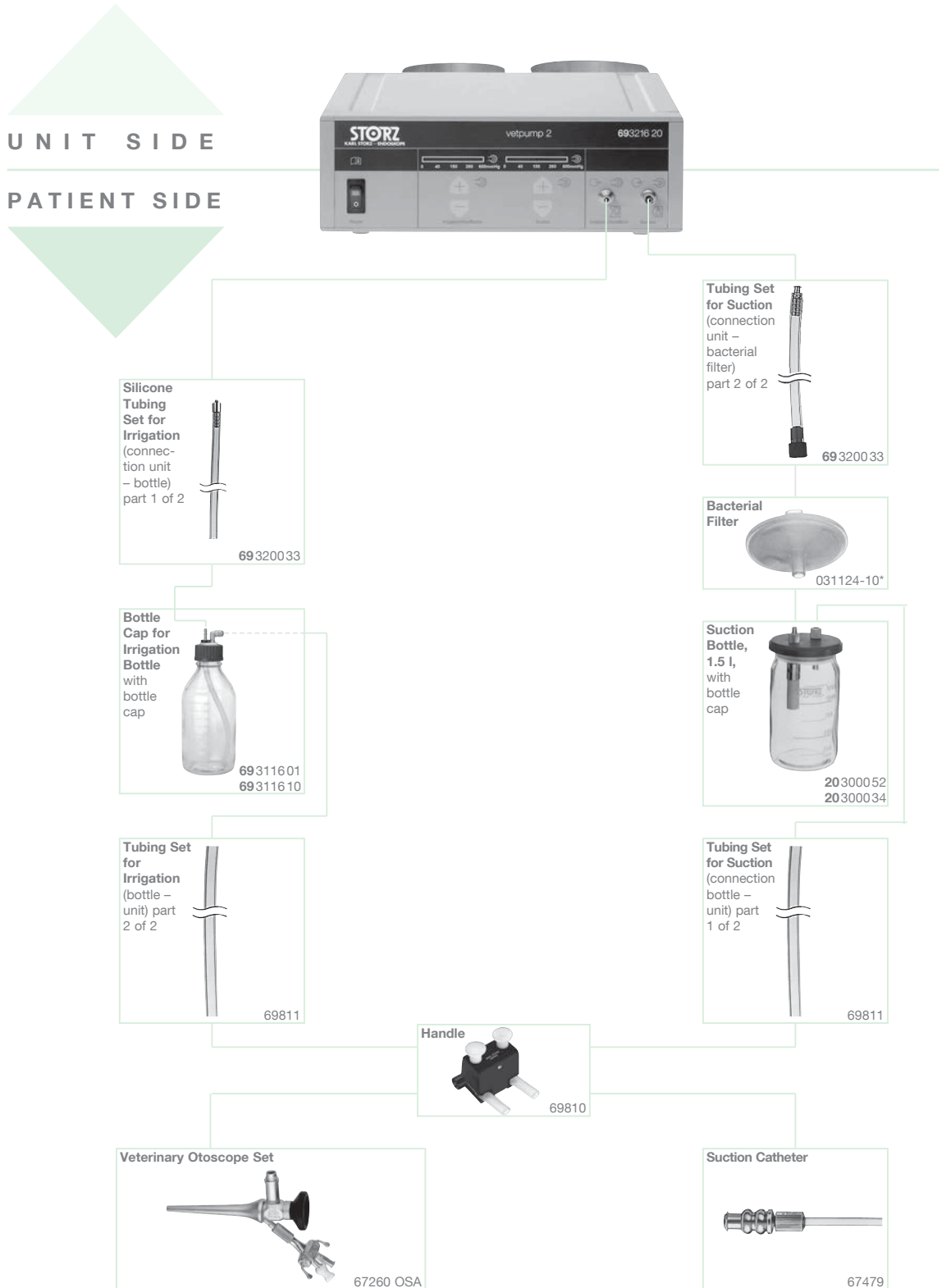


*



mtp medical technical promotion gmbh, take-off GewerbePark 46,
78579 Neuhausen ob Eck/Germany, Tel.: +49 (0)7467 94504-0, Fax: +49 (0)7467 94504-99,
E-Mail: info@mtp-tut.com, www.mtp-tut.com

VETPUMP® 2

System Components for use in Laparoscopy / Thoracoscopy

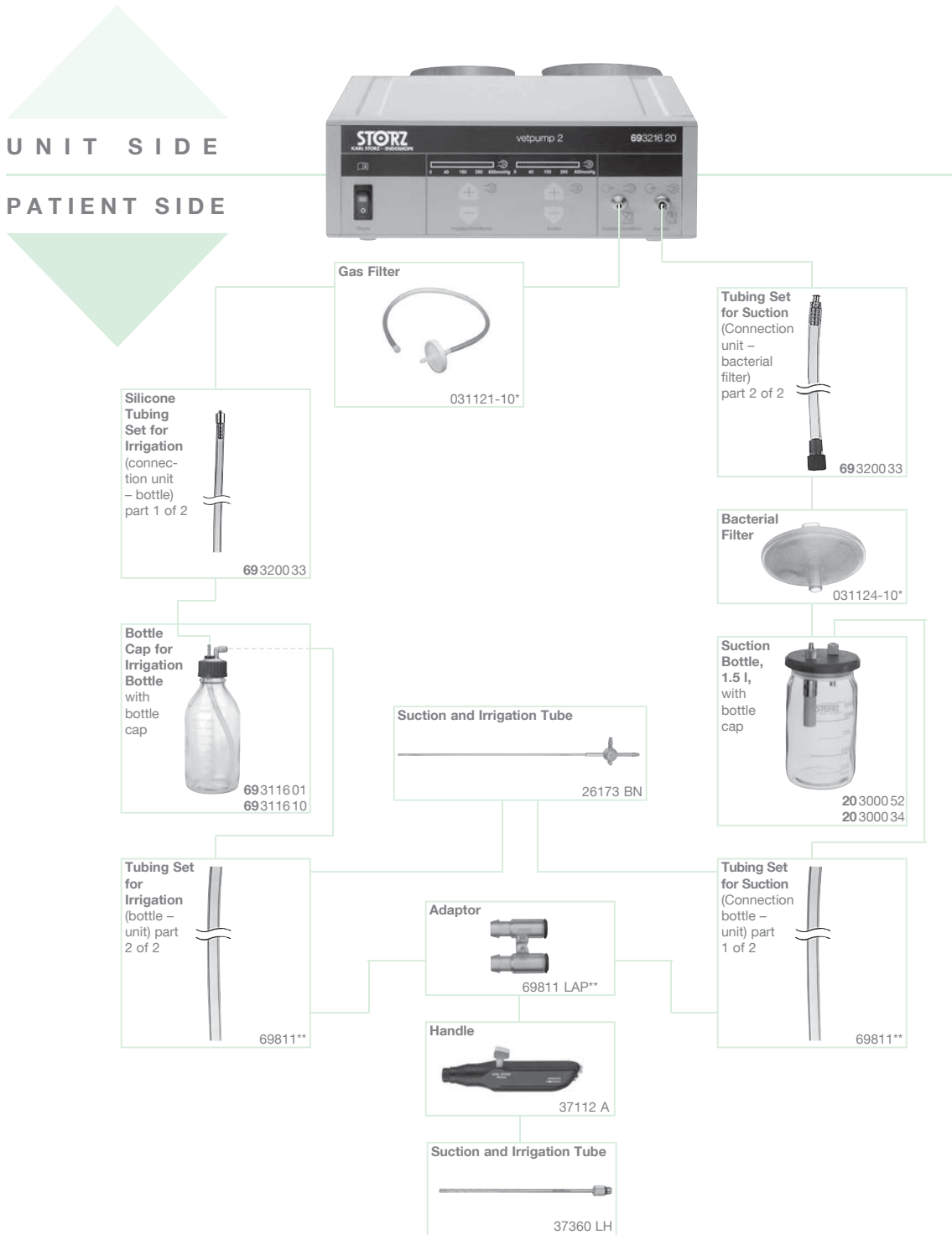

**Remove black plastic fitting from Tubing Set 69811 and replace with Adaptor 69811 LAP.


*



mtp medical technical promotion gmbh, take-off GewerbePark 46,
78579 Neuhausen ob Eck/Germany, Tel.: +49 (0)7467 94504-0, Fax: +49 (0)7467 94504-99,
E-Mail: info@mtp-tut.com, www.mtp-tut.com

VETPUMP® 2

System Components for use in Gastrointestinal Endoscopy

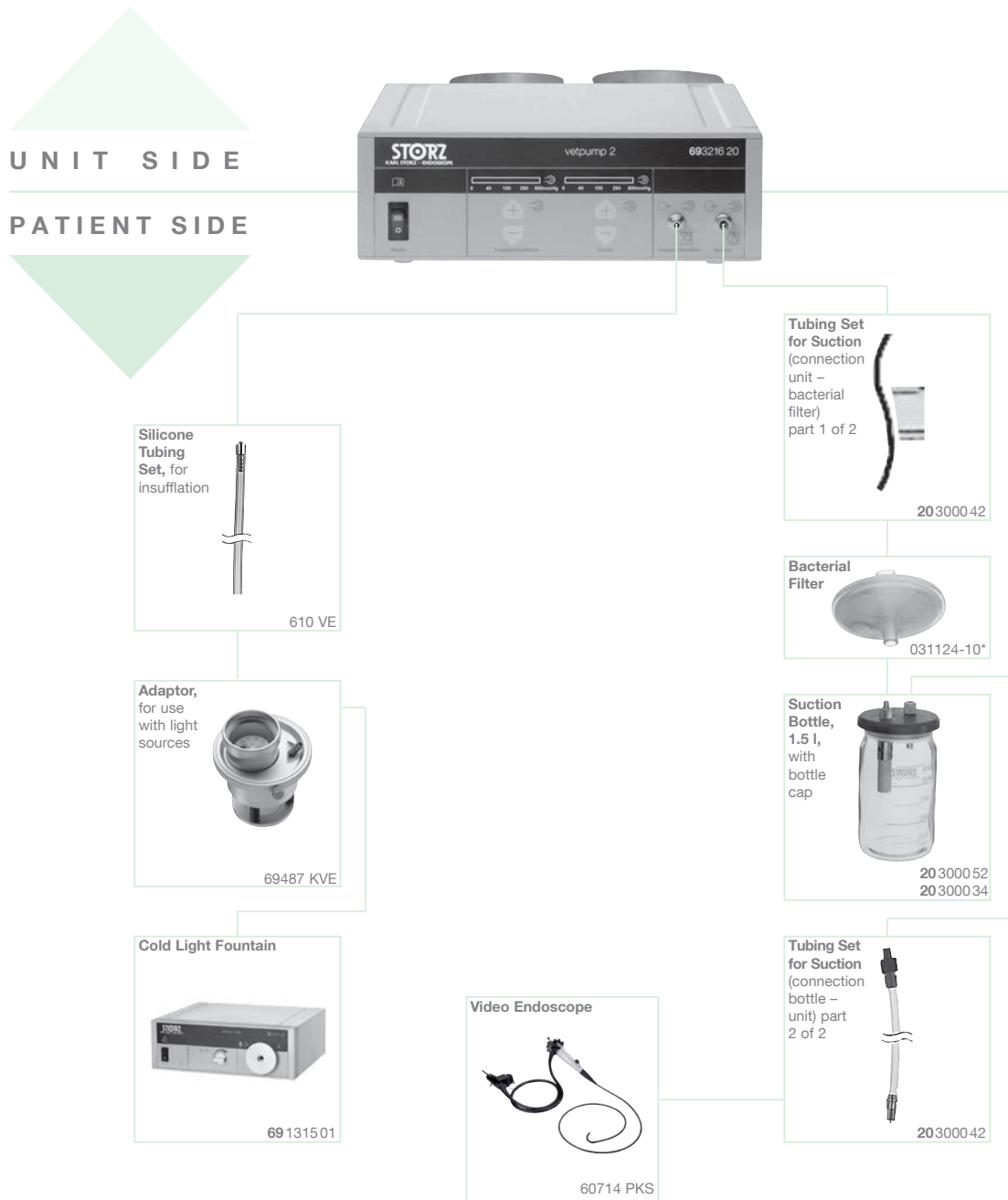

** Remove male LUER Lock fitting from end of Tubing **20300042** and connect free end to suction connector of endoscope.

*



mtp medical technical promotion gmbh, take-off GewerbePark 46,
78579 Neuhausen ob Eck/Germany, Tel.: +49 (0)7467 94504-0, Fax: +49 (0)7467 94504-99,
E-Mail: info@mtp-tut.com, www.mtp-tut.com



KARL STORZ GmbH & Co. KG
Mittelstraße 8, 78532 Tuttlingen, Germany
Postbox 230, 78503 Tuttlingen, Germany
Phone: +49 (0)7461 708-0
Fax: +49 (0)7461 708-105
E-Mail: vet@karlstorz.com

www.karlstorz.com